




SMART  
BEAUTY

2024-25

# THE BEST BEAUTY PERFOR MANCE.

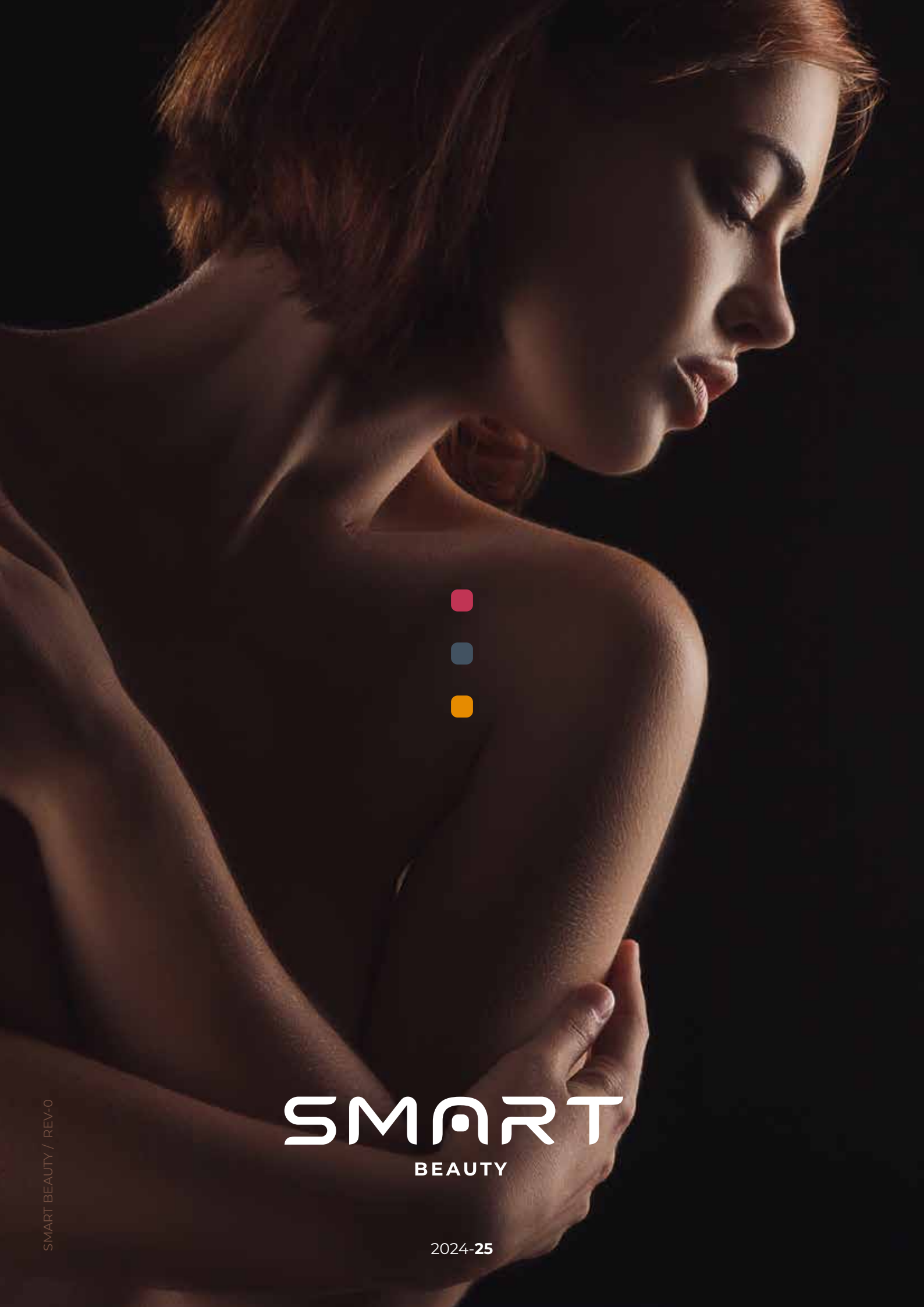

# SMART

BEAUTY



## ENJOY THE BEAUTY

Beauty has never been so accessible and performing as with the Smart Beauty new generation equipment.

SIMPLE   MULTIFUNCTIONAL  
ACCESSORIZED RESOLUTIVE  
TECHNOLOGICAL

Smart Beauty is your choice of advanced aesthetics. Made in Italy, high tech, compact, functional and efficient, it's the answer to any imperfection of face and body. It's the innovation you've been waiting for.

# FRT

FACE REPAIR TREATMENT



# HIGH TECHNOLOGY AT YOUR SERVICE

Smart Beauty is the global  
approach to facial beauty.

Smart Beauty offers you different multifunctional models to perform complete facial treatments as well as mono-function models to complete the equipment of your beauty center with a specific technology.


Very performing, designed and made in Italy, Smart Beauty devices for professional aesthetics are equipped with the most advanced technologies to give an effective response to the imperfections that cloud the splendor of the face and décolleté, with visible results in very few sessions.

These devices are very versatile and collaborative, the compactness of the dimensions, in fact, allows their placement even on small shelves, while an advanced software helps you to set up and manage the programs and steps of the treatments by viewing everything on an intuitive 7" monitor.

With Smart Beauty you can rely on devices complying with the provisions of the law, guaranteed 12 months, and certified ISO 9001, ISO 13485 and ISO 45001 by Bureau Veritas, to enable you to work with the utmost tranquility, and ensure safety to your customers.

# ANTI-AGE

## FACE REPAIR ANTI-AGE

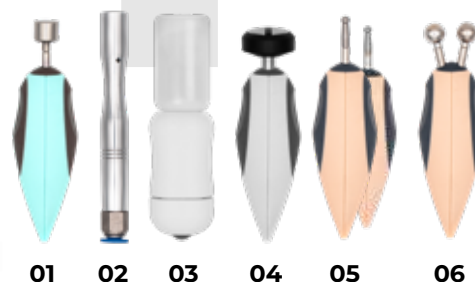

SMART.800-6-N  
dim. W.470 D.430 H.250 mm



Optional Cart  
dim. W.520 D.520 H.1120 mm

## Anti-age treatment

Time is not the only cause of skin aging, environmental factors such as pollution, and light in particular, can adversely affect the skin of the face, making it dry, dehydrated, not very elastic, with the classic and unattractive imperfections. Remove the impurities, deeply cleanse, nourish and tone the skin are important actions to brighten the face and preserve the beauty capital.

With Face Repair anti-age you can work on the face and décolleté of mature skin showing laxity and wrinkles, as well as on skin that's only at the beginning of the aging process. The result will be in any case exciting: toned, hydrated, oxygenated and lifted skin, smooth and compact texture, glowing face.



01



02



03

## 01 MOISTURE CONTROL

A non-hydrated skin undergoes premature aging. It is therefore important to check the level of hydration to evaluate how to intervene and rebalance it. A pre-treatment step that takes place with the Moisture Control function, a light and delicate current passage, that allows you to better determine the facial skin condition.

## 02 MICRODERMABRASION

The delicate and deep exfoliating action of Diamond Microdermabrasion, made with innovative diamond tips, effectively removes dead cells, toxins, impurities, promotes the cellular renewal of the treated tissue and the natural production of collagen. Improves blood circulation and skin permeability, facilitating greater absorption of the cosmetics chosen for treatment.

## 03 ELECTROPORATION

Electroporation is a skin tissue biorevitalization process. A transdermal conveying technique acting directly on cellular mechanisms, through modulable electrical impulses and vibrations generated by special dynamic currents. It is performed using a handpiece that exerts a weak electric discharge that frees the pores thus allowing the specific active ingredients to better penetrate, obtaining a result on the imperfection to be treated, favoring regeneration and rejuvenation. You'll see firmer, fuller, smoother skin from the first session.

## 04 RADIOFREQUENCY

Radiofrequency is effective to counteract skin relaxation, wrinkles and signs of aging, for a natural and progressive skin lifting effect. It acts through Diathermy, which reactivates fibroblasts for the production of hyaluronic acid, and collagen and elastin fibers, induces an increase in the amount of oxygen that is transported to the tissues, and speeds up cellular metabolism and the elimination of catabolic waste.

## 05 LIFTING

The Electro-lifting exploits the action of particular currents that simulate physiological nerve impulses to obtain the restoration of tissues and their revitalization. This technology exerts a stimulation at the level of skin microcirculation, increasing tissues nourishment and oxygenation. It also empowers the action of fibroblasts by increasing the quantity and quality of synthesized collagen, and stimulates the lymphatic system and its functions.

## 06 STIMOL

This treatment gently and selectively tones the atone muscles, counteracting skin relaxation. Using a particular handpiece with spherical emitters, pulses are transmitted that exploit a waveform and a specific frequency for the stimulation of the facial muscles, making the treatment pleasant. The handpiece will be placed on certain points for the correct micro-contraction, making the stimulation targeted, selective and effective.




**Optional:  
Multipolar Radiofrequency**



FACE REPAIR  
YOUNGER



## FACE REPAIR YOUNGER

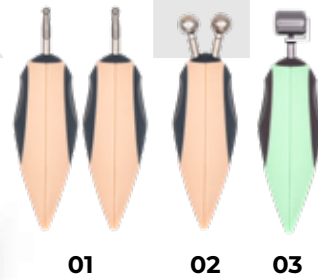

SMART.800-3  
dim. W.380 D.300 H.180 mm

## Facial rejuvenation specific treatment

Tired skin? Needs to regain radiance and youth? With Face Repair Younger, in a few steps, you can work on mature skin with visible signs of aging such as laxity and wrinkles.

In a short time, the skin will regain nourishment, hydration and oxygenation, will be toned, and with a smoother and more compact texture. In a word: lifted.

Simple and  
Intuitive  
Software



# RE

## 01 LIFTING

The Electro-lifting exploits the action of particular currents that simulate physiological nerve impulses to obtain the restoration of tissues and their revitalization. This technology exerts a stimulation at the level of skin microcirculation, increasing tissues nourishment and oxygenation. It also empowers the action of fibroblasts by increasing the quantity and quality of synthesized collagen, and stimulates the lymphatic system and its functions.

## 02 STIMOL

This treatment gently and selectively tones the atone muscles, counteracting skin relaxation. Using a particular handpiece with spherical emitters, pulses are transmitted that exploit a waveform and a specific frequency for the stimulation of the facial muscles, making the treatment pleasant. The handpiece will be placed on certain points for the correct micro-contraction, making the stimulation targeted, selective and effective.

## 03 IONTOPHORESIS

With Iontophoresis, the active ingredients used to take care of a specific imperfection are conveyed into the deep layers of tissues. A targeted and effective treatment, which contributes to the optimization of results in facial rejuvenation.



# SCRUB BRUSH



## FACE REPAIR SCRUB BRUSH DIAMOND




01

02

SMART.825  
dim. W.380 D.300 H.180 mm


Simple and  
Intuitive  
Software

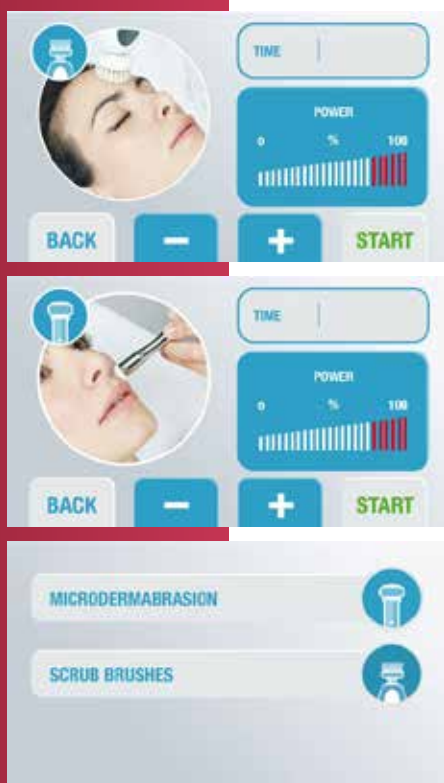

## Deep cleansing specific treatment


A good cleansing is essential to always have a young, healthy and radiant complexion. Deep cleansing and peeling treatment frees skin from toxins and impurities, and helps it to regenerate and “breathe” better as it stimulates microcirculation and the production of collagen, indispensable for its health. With Face Repair Scrub Brush Diamond you can with no problem work on all skin types, obtaining skin brightness and well-being outstanding results. Particularly useful as a pre-treatment, as it facilitates a better absorption of active ingredients, its use is particularly effective with oily and mixed skin, where it is necessary to rebalance the production of sebum.



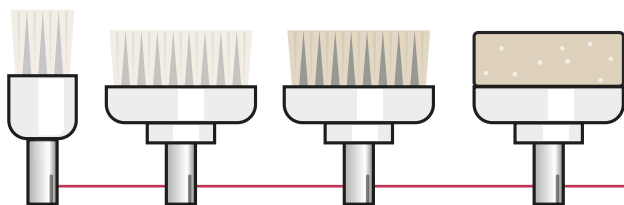

## 01 SCRUB BRUSHES

The Scrub Brushes treatment is a gentle mechanical peeling performed using different types of brushes, that can be operated at the desired speed. An ideal treatment for a deep cleansing of the skin. Its use is also recommended before the application of the cosmetic product, to facilitate the absorption of active ingredients.

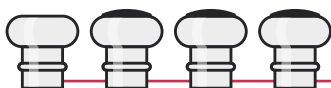

## 02 MICRODERMABRASION

The delicate and deep exfoliating action of Diamond Microdermabrasion, made with innovative diamond tips, effectively removes dead cells, toxins, impurities, promotes the cellular renewal of the treated tissue and the natural production of collagen. Improves blood circulation and skin permeability, facilitating greater absorption of the cosmetics chosen for treatment.



01



02



02

# ULTRA SCRUB BRUSH



SMART.814-3-N  
dim. W.380 D.300 H.180 mm

Simple and  
Intuitive  
Software



## FACE REPAIR ULTRA SCRUB BRUSH




### Deep cleansing specific treatment.

Whatever the age and skin type, facial cleansing is an important treatment to take care of the beauty capital of both women and men, to maintain healthy skin when you are young and to prevent skin aging in adulthood. The thorough cleansing of the face that you can do with Face Repair Ultra Scrub Brush, in addition to freeing the skin from impurities, toxins and dead cells, acts on the microcirculation and reactivates it, bringing more oxygenation and nourishment. In addition, it immediately gives a bright complexion and allows the skin to better absorb cosmetic products to treat imperfections, giving a pleasant "breathing skin" sensation.




## 01 25kHz ULTRASOUND

Deep cleansing action using a high frequency vibrating spatula to perform an effective removal of surface dead cells. The vibrations reactivate the microcirculation and stimulate the production of collagen, elastin, and endogenous hyaluronic acid. The use of different frequencies and 14 specific programs allows you to perform a regenerating, toning and elasticizing action, with excellent results on all skin types.

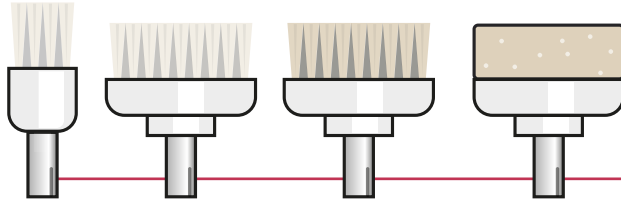

## 02 SCRUB BRUSHES

The Scrub Brushes treatment is a gentle mechanical peeling performed using different types of brushes, that can be operated at the desired speed. An ideal treatment for a deep cleansing of the skin. Its use is also recommended before the application of the cosmetic product, to facilitate the absorption of active ingredients.

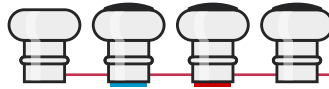

## 03 MICRODERMABRASION

The delicate and deep exfoliating action of Diamond Microdermabrasion, made with innovative diamond tips, effectively removes dead cells, toxins, impurities, promotes the cellular renewal of the treated tissue and the natural production of collagen. Improves blood circulation and skin permeability, facilitating greater absorption of the cosmetics chosen for treatment.



# MICRODERMABRASION

## FACE REPAIR MICRODERMABRASION



SMART.820-N  
dim. W.380 D.300 H.180 mm


Simple and  
Intuitive  
Software



### Cell regeneration specific treatment.

The Microdermabrasion technique is a treatment that promotes cell regeneration. A non-invasive method, which aims to a smoother skin of the face, reducing wrinkles, signs of aging, scars and other imperfections.

With Face Repair Microdermoabrasion you can perform a more soft and skin-friendly abrasive action. It is an ideal solution as a pre-facial treatment, but it is also very useful when the skin lacks brightness or when it is afflicted by more important imperfections such as couperose and skin spots.

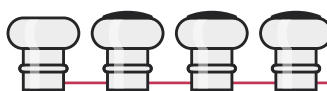

### MICRODERMABRASION

The delicate and deep exfoliating action of Diamond Microdermabrasion, made with innovative diamond tips, effectively removes dead cells, toxins, impurities, promotes the cellular renewal of the treated tissue and the natural production of collagen. Improves blood circulation and skin permeability, facilitating greater absorption of the cosmetics chosen for treatment.

# ULTRASOUND

## FACE REPAIR ULTRASOUND



SMART.814-N  
dim. W.380 D.300 H.180 mm

Simple  
and Intuitive  
Software



## Deep cleansing specific treatment

Cleansing, exfoliating, hydrating; facial cleansing is a very important basic treatment to prevent or treat skin imperfections. Face Repair Ultrasound acts on the skin surface through a fast vibratory movement, triggering the exfoliation of skin surface cells as well as the removal of sebum and impurities of environmental origin, such as fine dust and smog. It reactivates cell turnover and promotes the action of fibroblasts responsible for the production of collagen. The result is a radiant and deeply purified skin, an immediate "flash" effect. Ideal for thick skin that need deeper cleansing, the facial cleansing with ultrasounds is not suitable in the presence of skin inflammatory pathologies, neoformations, atopic dermatitis or rosacea.

### 25kHz ULTRASOUND

Deep cleansing action using a high frequency vibrating spatula to perform an effective removal of surface dead cells. The vibrations reactivate the microcirculation and stimulate the production of collagen, elastin, and endogenous hyaluronic acid. The use of different frequencies and 14 specific programs allows you to perform a regenerating, toning and elasticizing action, with excellent results on all skin types.



# ELECTROPORATION



SMART.755-N

dim. W.380 D.300 H.180 mm

## FACE REPAIR ELECTROPORATION



### Specific treatment to enhance active ingredients action.

You identified the imperfection to be treated. You analyzed the quality of the skin. You already decided which active ingredient to use to cope with imperfection. Now you need to proceed to treatment, but how can you speed up the beauty effect you want to obtain?

With Face Repair Electroporation, you can make the active ingredients better penetrate into the deeper layers of the skin, and more effectively act on the imperfection.

#### ELECTROPORATION

Electroporation is a skin tissue biorevitalization process. A transdermal conveying technique acting directly on cellular mechanisms, through modulable electrical impulses and vibrations generated by special dynamic currents.

It is performed using a handpiece that exerts a weak electric discharge that frees the pores thus allowing the specific active ingredients to better penetrate, obtaining a result on the imperfection to be treated, favoring regeneration and rejuvenation. You'll see firmer, fuller, smoother skin from the first session.



# STIMOL

## FACE REPAIR REGENERATION



SMART.824-N  
dim. W.380 D.300 H.180 mm

Simple and  
Intuitive  
Software



## Specific Stimulating treatment

Do you need to give a boost of beauty to face and décolleté skin? With Face Repair you can offer your customers a quick and simple solution to give new vitality and youth to the skin. A treatment that can deeply act, for a lifting effect visible from the first sessions.

### BIPOLAR RADIOFREQUENCY

This treatment gently and selectively tones the atone muscles, counteracting skin relaxation. By raising the temperature of the existing connective tissue and, through heat, activating the fibroblasts to produce new collagen and new elastin, Bipolar Radiofrequency allows to mitigate wrinkles and signs of aging, for a whole natural compacting and smoothing effect on tissues.

**Optional:**  
**Multipolar Radiofrequency**



# VAPORIZER

020

## NEW PHYTO VAP



## VAPORIZER

Need to work on acne, comedones, asphyxiated skin, impurities? With the New Phyto Vap vaporizer you can open pores, eliminate toxins, clean skin, moisturize and improve blood circulation. And not only that.

Thanks to the three different baskets in which you can insert, and change when required, the essential oils useful for the treatment you are performing, New Phito Vap can optimize its effectiveness, as well as make the beauty session pleasant and beneficial.



# VAPORIZER

## NEW OVER VAP



Included:  
code 9385 Boiler  
code 9004-N2 Stand

Included:  
code 9385 Boiler  
code 9004-N2 Stand

The New Over Vap vaporizer has been engineered to match both ozone and aromatic herbs or essential oils with the produced water vapor. Ozone is known for its antibacterial action and its ability to affect the acidity of the skin. Aromatic and phytocosmetic herbs are ideal for treating specific imperfections.

In the steam jet, thanks to the tool on the outlet nozzle, these essential substances slowly release their active ingredients, spreading on the epidermis in a microscopic way. This allows you to optimize their beneficial functions.

code 9385  
Boiler



code 9395  
Herb baskets

code 9390 Essence  
holder spare part



# EA

EPIL ACTIVE TREATMENT



## TOP TECHNOLOGY

A complete range of products for unwanted hair removal.

Hair removal has become an essential part of body care routine for women and men. One of the most requested treatments in beauty centers.

Smart Beauty gives you the opportunity to choose between different hair removal systems, to respond to any customer request and expand the range of services in your beauty center.

Laser Diode devices to enrich your equipment with the most advanced and best technology for permanent removal of unwanted hair, flash needle to work in particular areas of the body and for finishing, and then wax warmers to perform waxing, one of the most popular hair removal methods.

Our products are manufactured strictly respecting the regulations that guide our industry, according to the information sheets of Decree no. 206/2015, and guaranteed 12 months, thus allowing you to run your business in total serenity and safety, giving tranquility to those who rely on you to solve the problem of unwanted hair.

LASER

EPILATION ACTIVE  
LASER



2188-06-64-SMA  
dim. W.470 D.430 H.250 mm

## 808nm Laser Active Hair Removal Treatment

A velvety skin without unwanted hair? You can give a dream to your customers, with the Laser active hair removal by Smart Beauty. Targeted, fast, painless, and comfortable, with the advanced Laser Diode Technology you can work with a safe method on all parts of the body and treat any type of hair, even the harder and more difficult to eliminate male ones. Thanks to its wavelength (808nm), the light beam emitted by the laser warms the hair stem and its root, weakening it. In this way, session after session, there is a progressive thinning and a significant slowing down of the regrowth, until the disappearance of the hair. The effectiveness of the system is proven by the numerous tests carried out by the company.





dimensions W.520 D.520 H.1120 mm

## 5 SPOTS PER SECOND

The Laser 808nm system of Smart Beauty allows you to greatly improve the treatment times and effectiveness. With the multi-spot methodology, through its ergonomic handpiece, equipped with an anti-fall strap, laser emission lights and cooling lights, you can dispense up to 5 spots in the fraction of a second. A great advantage for you who want safety, precision of work, results, and speed in execution.

## HD MATRIX

You will experience the great innovation of this hair removal system: the 600 Watts HD-MATRIX technology. This technology allows to hit the hair bulb in a more targeted way because the emission of the high-definition laser beam is much more uniform and does not produce dispersion. The Laser wave, absorbed by the melanin of the hair stem, turns into heat, increases the temperature, and triggers a natural physiological process that leads to the progressive disappearance of unwanted hair.

## MAINTENANCE

Facing such advanced equipment you are surely wondering how much its maintenance can affect. Well, this hair removal system will surprise you because it requires minimum maintenance every 36 months. In the instruction booklet you will find instructions for daily cleaning and maintenance of the equipment.



MAX.40J/cm<sup>2</sup>-808 nm 300ms MAX EN 60825-1



# NEEDLE

## EPILATION ACTIVE NEEDLE




SMART.501-N  
dim. W.380 D.300 H.180 mm

## Flash Progressive Hair Removal

The hair removal system with the flash needle offers you the possibility to work on unwanted hair in a gentle, progressive and painless way, until they are eliminated. Using the Thermolysis technique, a small needle inserted under the skin hits the hair bulb, destroying it through the use of electricity. It can be used on larger and smaller areas of the body and is effective on hairiness, so it is also ideal as a finish, when unwanted hair, treated with the Laser, become thinner.

**Pulse emission times from 0.125 to 2 seconds**

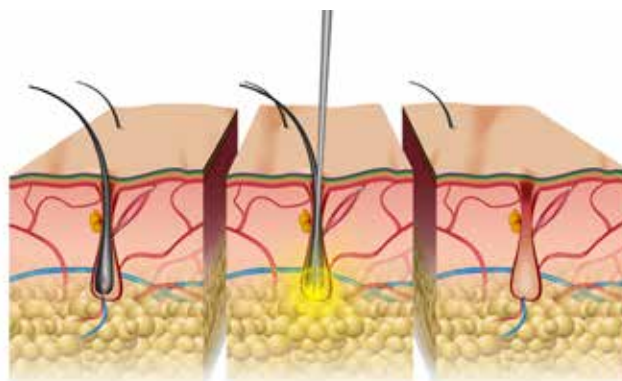



## ISOLATED



**F3.9191**  
0,075 mm  
facial hair



**F4.9192**  
0,102 mm  
medium hair



**F5.9193**  
0,127 mm  
medium hard hair

## NOT ISOLATED



**F3.9177**  
0,075 mm  
facial hair



**F4.9178**  
0,102 mm  
medium hair




**F5.9179**  
0,127 mm  
medium hard hair

## 6 DIFFERENT NEEDLES

The needles, made of stainless steel of the highest quality and in various sizes, allow you to work on any type of hair ensuring the constant supply of power for the duration of treatment. The diamond-drawn tip allows the perfect insertion, while the stem is isolated by a plastic film for maximum safety.

## ADVANCED SOFTWARE

The system is managed by an advanced software, you can easily adjust the power and duration of thermolysis and see it all on the 7" display.

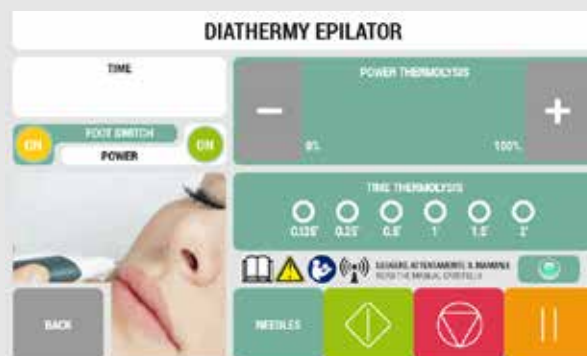

# WAX HEATERS



**OVER2200N-500** ●  
(BLACK 400-500 grams)

**OVER2200N-800** ●  
(BLACK 800 grams)

**OVER2210B-500** ●  
(WHITE 400-500 grams)

**OVER2210B-800** ●  
(WHITE 800 grams)



## LUNA - MONO - TRIO

### Wax heaters

Hot waxing is among the last minute methods to have a smooth skin and free from unwanted hair . A treatment that is not definitive but allows your customers to be neat for at least four weeks. This is why a tool that can not miss in your salon is the wax heater. Smart Beauty offers you three models designed to meet your needs and those of your customers, but be careful not to treat those suffering from capillary fragility or varicose veins.

### LUNA

Compact, practical and functional, the Smart Beauty Luna wax heater, suitable for all types of wax, is available in two versions, one with a 400 grams jar holder and the other with a 800 grams jar holder. The thermostat allows you to easily adjust the temperature from 30° C to 120° C, to reach the desired temperature and keep it constant. With a sinuous design, you can choose between two colours, white and black, to better harmonize inside your beauty centre.

**EC DECLARATION OF CONFORMITY**  
**FOR ITALY, IT IS HEREBY DECLARED THAT THE**  
**CHARACTERISTICS REQUIRED BY THE TECHNICAL**  
**INFORMATION SHEET No. 10 ARE RESPECTED.**

**List of Devices: Wax Heater for Waxing**  
**DECREE no. 206 of 15 October 2015.**



210



215



218



219



## TRIO

The use of honey waxes allows you to make gentle depilatory treatments in any part of the body. For this, the wax should be heated and kept at a constant temperature throughout the duration of treatment. Smart Beauty Trio wax warmer is what you were looking for. With the thermostat you can adjust the temperature from 30°C to 120°C for the wax jar, while the cartridge thermostat is adjustable from 20°C to 85°C. Reliable, practical and safe, it includes a universal cartridge holder and the insertion of a 400g/500g and 800g wax jar. Also available with mono wax compartment.

## MONO

Small yet very performing, Smart Beauty Mono Soft Touch wax heater has a double resistance to quickly reach the desired temperature and is equipped with a double window to display the cartridge fill level. It's built with nylon plastic outer materials, easily cleaned with oil, and has an ergonomic handle that facilitates manual dexterity, designed to prevent overheating of the hand. Also proposed in the 218 version that allows you to adjust the temperature and then to work with all types of wax.

### TECHNICAL FEATURES

1. Double resistance to quickly reach the desired temperature
2. Nylon plastic outer materials (easily cleaned with oil, without damaging the wax heater)
3. Double window to see cartridge fill level
4. Ergonomic handle that facilitates manual handling, and designed to be kept at a temperature that prevents operator's hand overheating
5. The wax heater shall be mounted under pressure and with the aid of four fixing screws, to prevent breakage due to misuse
6. Includes power cord inside
7. The temperature adjustment is available in the 218 version, this feature makes this object unique and suitable for all types of wax.

Dimensions in mm: W. 350 D. 300 H. 160  
120/240V · 20W (120W max) Class II · 50/60Hz.



# BS1

BODY SLIM TREATMENT

# INNOVATIVE SYSTEMS FOR ADVANCED AESTHETICS

The best answer for body  
beauty without surgery.

When the goal is to keep a toned physique, intervene on critical points by reducing fat and orange peel, and reshaping and slimming the silhouette, the answer is Smart Beauty. A beauty that you can conquer without a scalpel, through personalized, non-invasive treatments, which do not provide recovery time, thanks to the high technology of the latest generation of our equipment.

Safe, certified, 12-month warranted devices. Technologies that comply with the regulations in accordance with the technical information sheets attached to Decree No. 206/2015.

A wide range of solutions to complete and enrich the endowment of your Institute. Simple use and high performance devices, to perform professional body treatments with the certainty of always offering your customers the best service.

Compact and small in size, Smart Beauty devices can fit in any room, allowing you to optimize space.

# SWINGING MASSAGE



SMART.3190

dim. W.470 D.430 H.250 mm



Optional Cart

dim. W.520 D.520 H.1120 mm

## BODY SLIM SWINGING MASSAGE



### Continuous alternating pressure Swinging Massage

Rolls, weighted shapes, sagging skin, cellulite\*, fat\* and atony, with Smart Body Slim Swinging Massage you can solve these problems with a very effective continuous alternating pressure Shock Wave treatment. Also perfect for customers who ask for a last minute treatment to get in shape.

The results? Reduction of cellulite\* and fat\*, reshaped body, toned and hydrated skin, and a more compact texture, already from the first sessions.

- SIMPLE TO USE
- REDUCES CELLULITE\* IMPERFECTION
- IMPROVES SKIN TONE AND ELASTICITY
- STIMULATES COLLAGEN PRODUCTION
- VISIBLE RESULTS IN A SHORT TIME
- REDUCES FAT\* PANNICLES
- RESULTS ALSO ON LOCALIZED FAT\*
- STIMULATES BIOLOGICAL PROCESSES
- NO RECOVERY TIME

\*Cellulite: imperfections caused by cellulite

\*Fat: adjuvant in the treatment of localized fat

## SHOCK WAVE


The Shock Wave is a real biostimulation treatment that reactivates the normal physiological processes of tissue repair and regeneration, increases veno-lymphatic microcirculation, reactivates metabolism, encourages lipase, and tones. Without thermal effects, without pain, quickly and individually, you can act in maximum safety with a strong wave of alternating pressure that, crossing the fluids, acts only and exclusively on adipocytes, remodeling fat deposits and cellulite.

## HIGH TECHNOLOGY AND DESIGN

The equipment is managed by an advanced software that guarantees very advanced solutions and an extremely intuitive interface. The 7" display facilitates easy access to controls, programming, and immediate work checking. The compact design makes this equipment easily positioned in all working spaces, even in small cabins.

## SWINGING MASSAGE

### EFFECT

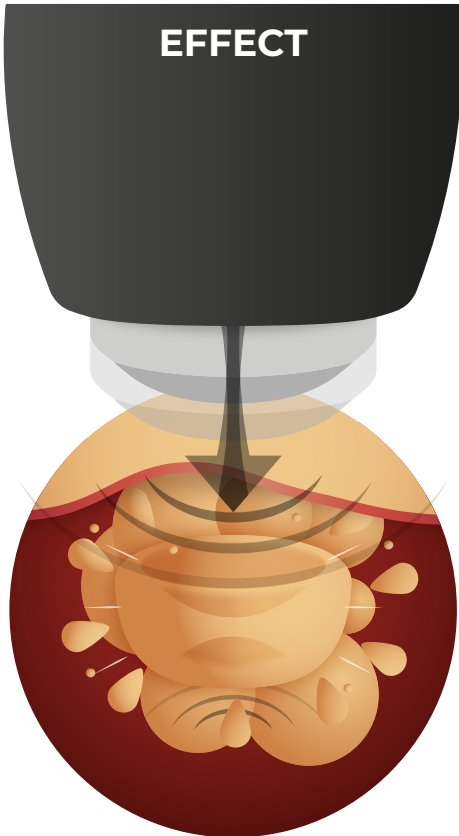

# INFRARED

## BODY SLIM INFRARED



SMART.719-N  
dim. W.470 D.430 H.250 mm



Optional Cart  
dim. W.520 D.520 H.1120 mm

Simple and  
Intuitive  
Software



## Infrared

One of the best methods of heat transmission is infrared radiation. The treatment pleasantly warms the part, penetrating into the adipose tissue, the localized increase in temperature promotes vasodilation and therefore a correct veno-lymphatic flow that facilitates the drainage of fluids and toxins. Thanks to a particular wavelength, very selective, the photobiological radiation emits the heat only when it reaches the adipose panniculus, increasing the temperature and the enzymatic activity consequently, favoring local lipolysis, splitting the storage triglycerides into fatty acids and glycerol to enter them in the blood sector and be metabolized and consumed as an energy source by contractile muscle activity. INFRARED SMART is a computerized electronic equipment programmed for thermo-localized treatments for aesthetic applications, aimed at achieving the ideal silhouette.


The heart of the INFRARED SMART is a computerized central unit entirely electronic, equipped with a special microprocessor able to organize, following the indications given by the aesthetic operator, the thermal energy applied to the areas affected by the imperfection of fat\* or cellulite\*, favoring the obtaining of slender and shaped thighs, toned buttocks and abdomen, and a thinned waist.

\*Cellulite: imperfections caused by cellulite

\*Fat: adjuvant in the treatment of localized fat

The use of INFRARED SMART is very simple: the areas affected by imperfections are wrapped by special bands, the SLIMBANDS, specific for each sector of the body and packed with a special non-toxic material connected to the equipment, which guarantees its operation with a maximum tension of 24 Volts. Slimbands allow the flow of thermal energy that develops a pleasant and uniform heat, deep and localized, wisely modulated and controlled by the micro-processor that also guarantees absolute safety for the customer.

After a series of treatments with INFRARED SMART, you will notice a more elastic skin and improved epidermis. Even more noticeable will be the results in the reduction of fat\* deposits and cellulite\*, in addition to a decrease in weight and centimeters.

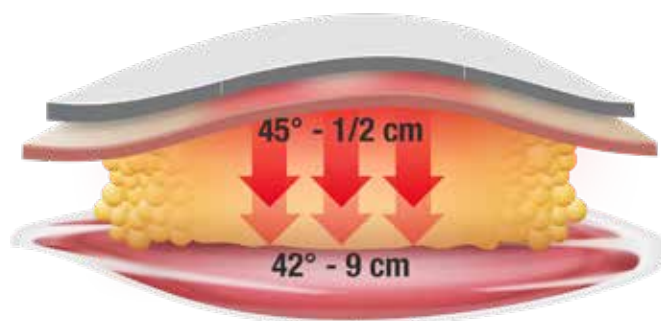

\*Cellulite: imperfections caused by cellulite  
\*Fat: adjuvant in the treatment of localized fat



# ENDOMASSAGE

## BODY SLIM ENDOMASSAGE



3238-34-31X-SMR  
dim. W.470 D.430 H.250 mm




Optional Cart  
dim. W.520 D.520 H.1120 mm

## Endodermal massage

Cellulite\*, stretch marks, water retention, are imperfections often caused by sedentary behaviour and an unbalanced lifestyle. Regain a toned and fit body, you can, acting with Endomassage directly on the compromised blood and lymphatic circulation.



\*Cellulite: imperfections caused by cellulite

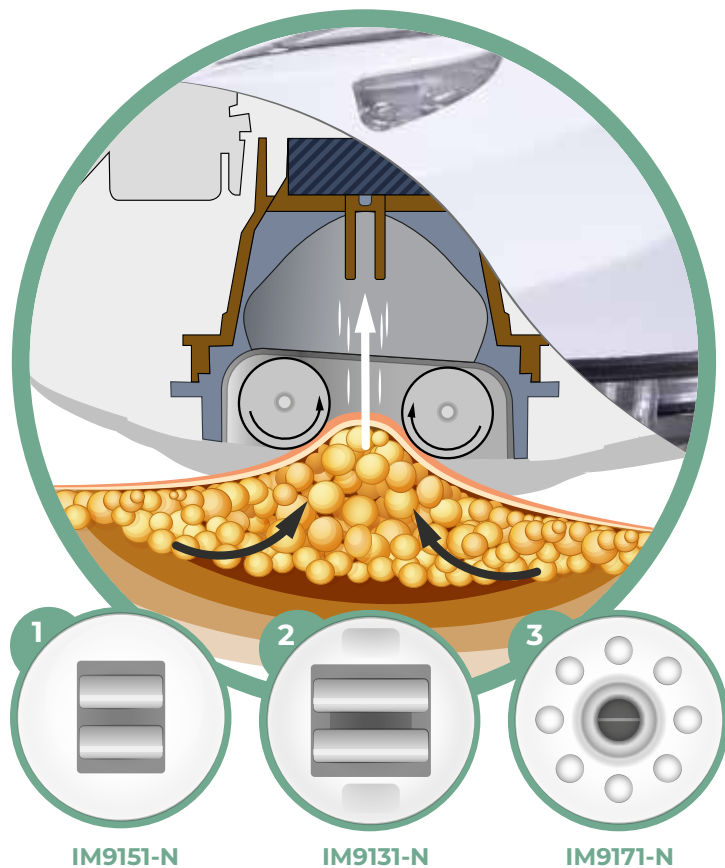

## ENDODERMAL MASSAGE

The best revitalizing treatment is the Endodermal Massage, a real aspirated connective massage, aimed at reactivating blood circulation and re-oxygenating the cells.

The special Vacuum-massage handpiece allows you to perform, simultaneously or individually, detachment, drainage and rhythmic massage treatments, with both automatic and manual working modes.

The handpiece uses the action of interchangeable rollers or spheres, which exert a positive pressure associated with the application of negative pressure on the skin and subcutaneous.

- 1 Accessory with 35mm rollers: for specific zones**
- 2 Accessory with rollers 55mm: for extended areas**
- 3 Accessory with spheres: to follow the shapes of the silhouette and make a reshaping massage**

Inside the handpiece there are special cleaning filters, which ensure maximum hygiene and longevity of the equipment. All accessories can be simply removed and replaced.



REDUCING

BODY SLIM  
REDUCING



752

dim. W.470 D.430 H.250 mm



Optional Cart

dim. W.520 D.520 H.1120 mm

## Slimmer with the 980nm Laser


Those extra inches that, despite the good habits, never go away. A problem for many people who attend your beauty salon.

So what to do? What to propose? To rid the body of excess fat, Smart Beauty presents Body Slim Reducing, a very effective and non-invasive slimming and remodeling equipment. A system that uses Laser LLLT 980nm technology for slimming treatment in association with a combination of special energies capable of acting effectively on the most common imperfections that undermine the beauty and harmony of the silhouette.



**LOW LEVEL LASER TREATMENT**  
**10 LASER EMITTERS INTEGRATED IN THE 980nm LASER**  
**PLATE**  
**STIMULATING SYNERGISTIC EMISSION**



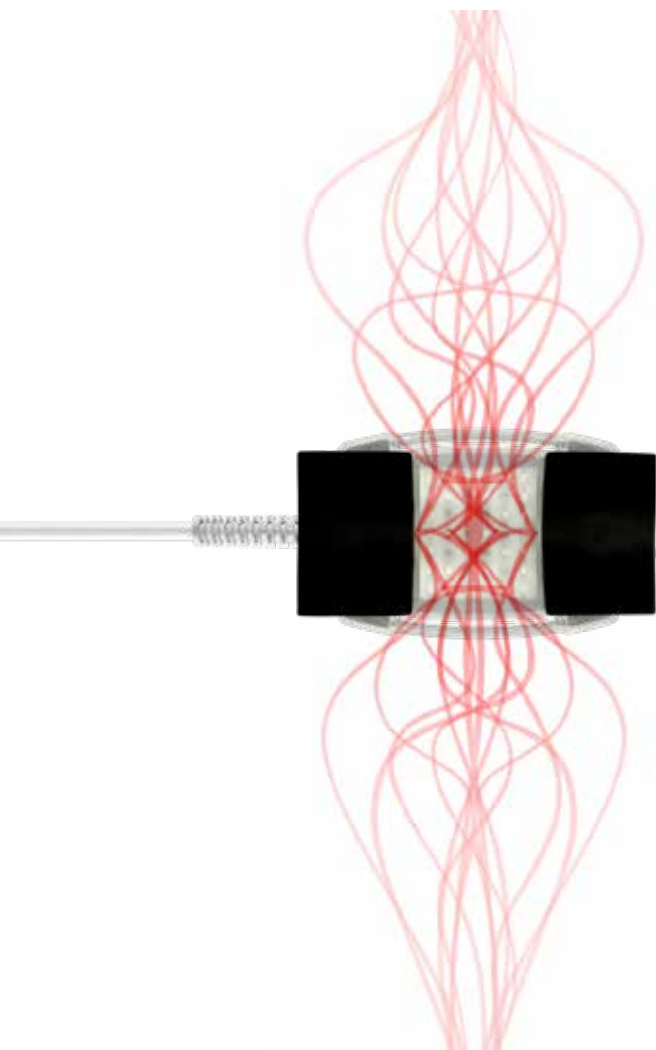

### **LASER LLLT**

An absolutely safe and latest generation treatment with the use of LLLT technology (Low Level Laser Technology) 980nm, suitable for those who want to get rid of excess fat\* without using needles and scalpel. You can treat different areas of the body, especially those where lipids settle, such as thighs, waist, hips, back, arms, abdomen.

The energy emitted by the laser affects the surface layer of fat\*, triggering the process of fat breakdown, which is then transported by water through the lymphatic system and kidneys, and then naturally expelled from the body.

### **8 PLATES - 16 PROGRAMMES**

The equipment, which works with 8 special plates that allow to treat several areas of the body simultaneously, allows to intervene on the different imperfections synergistically, with 4 paths divided into 16 different programs, designed ad hoc to obtain immediately visible results.

**Treatments are non-invasive and do not require recovery time before going back to your daily schedule.**

\*Fat: adjuvant in the treatment of localized fat



# RADIOFREQUENCY

## BODY SLIM RADIOFREQUENCY



SMART.745-N

dim. W.380 D.300 H.180 mm

OPTIONAL

## Radiofrequency

It is one of the most popular treatments for toning and remodeling the shapes, representative of a philosophy of regenerative and not transformative beauty, without the use of needles and scalpel; we are talking about Radiofrequency.

A technology that promises to give an increasingly compact, smooth and toned skin. Body Slim Radiofrequency is the equipment that can not miss in your Beauty Institute, to complete the offer of treatments.



# REF

OPTIONAL

- **Body Multipolar Radiofrequency**, exploits the characteristics of monopolar and bipolar methods to bring energy under the skin, focusing the range of action in a single point through the 4-pole system. In this way, higher temperatures can be obtained both on the surface and in the deeper layers.

- **Face Bipolar Radiofrequency**, a treatment with an extraordinary lifting effect that counteracts skin relaxation. By raising the temperature of the existing connective tissue and, through heat, causing the fibroblasts to produce new collagen and new elastin, it allows wrinkles softening, for a tissue compacting and smoothing completely natural effect.

- **Face Multipolar Radiofrequency**, acting through the increase of the subcutaneous temperature, which causes the denaturation of collagen proteins giving rise to a reorganization of the lines of skin tension. Moreover, it activates the fibroblasts that, as a result of the thermal insult, increase over time the production of collagen and elastin with a consequent increase in the consistency of the dermis, for a significant anti-age effect.

## BIPOLAR RADIOFREQUENCY

The bio-physical effect of Radiofrequency is based on the transformation of energy into heat. It increases local metabolism, speeds up blood flow in the microcirculation, eliminates fluid stagnation in the capillaries, frees interstitial tissue from edematous congestion. It also restores the necessary temperature to revive lipolysis, stops the fibrotic reaction of the subcutaneous tissue, reduces loose skin, and the effects of imperfections caused by cellulite\*.

The treatment is ideal to act on the imperfections of cellulite\*, fat\*, and skin laxity. Radiofrequency regenerates and revitalizes tissues in a painless and non-invasive way, giving the skin elasticity and tone, and can be applied on all parts of the body.

\*Cellulite: imperfections caused by cellulite

\*Fat: adjuvant in the treatment of localized fat



# PRESSOMASSAGE

## BODY SLIM PRESSOMASSAGE



SMART.732  
dim. W.470 D.430 H.250 mm



Optional Cart  
dim. W.520 D.520 H.1120 mm

## Presso-Massage

Do you want to help your customers to get rid of swelling, muscle tension, a few of extra pounds, and give them lightness, body harmony and relaxation for the mind? Smart Beauty presents three presso-massage equipment models, offering different programs to perform different types of peristaltic massages for legs and arms.

## The Cuffs

The cuffs are the flagship of this equipment. They provide an adjustable triangle to improve pressure at the level of the inguinal ganglia, and to manage the size at the height of the hips, they are also adaptable both in height and circumference, to ensure maximum versatility.

The progressive inflation and the internal and overlapping air chambers stimulate the veno-lymphatic return thrust, improve drainage and guarantee a gentle, deep and at the same time uniform massage on the treated parts, wrapping the limbs like a wave, from the bottom up.



**Code SMART.730**  
**Cart not available**



### **SMART PRESSO**

Smart Presso is the basic equipment that allows you to manage 3 programs for the legs and 1 for the arms. Improves the functioning of the circulatory and lymphatic system and is effective against imperfections caused by cellulite\*, and water retention.  
(Cart not available).

**Code SMART.732-W**  
**Optional Cart**



### **SMART PRESSO WARM**

Smart Presso Warm offers you 4 Presso-massage programs, 3 for the legs and 1 for the arms. This equipment combines the action of heat with these programs, thus enhancing the presso-massage treatment making it even more effective, and it also generates a reducing and at the same time pleasant action, a real pampering for the body.  
(Optional Cart)

**Code SMART.732**  
**Optional Cart**



### **SMART PRESSO SUPER WAVE**

Smart Presso Super Wave includes 6 Presso-massage programs, 4 dedicated to the legs and 2 to the arms. The device is equipped with a 8" display and specially designed to counteract the imperfections caused by cellulite\*, swelling, and water retention. Restores venous and lymphatic circulatory function by reducing swelling and muscle tension, while simultaneously slimming the silhouette and relieving stress.  
(Optional Cart).

\*Cellulite: imperfections caused by cellulite

### **Interchangeable Sectors**





# ACCESORIES SORIES

## SMART LENTE LED

Smart Lente Led is the ideal tool to better focus on imperfections or details that you want to see more clearly. Protected by a light fringe screen, with a biconvex magnifier in 3,5 diopter optical glass, this LED lens with cold light has a folding arm to adjust it according to needs, and an attachment for connection to the work table.

Optional: Stand with castors

**Code OVER9000-N**

SMART BEAUTY • 44



9001-N

9004-N1  
660 mm diam. stand, black,  
castors provided



# AUTOCLAVE 9910-N1

## TECHNICAL FEATURES

- Power supply 220 V, 1400 Watts, 50/60 Hz
- Class B
- Sterilization grades 121° - 134°
- Minimum sterilization time 4 minutes
- Maximum sterilization time 30 minutes
- Internal chamber dimensions 170 x 320 mm
- Chamber volume 9 Liters
- Dimensions 370 x 345 x 565 mm
- Weight 34,5 kg

The sterilization of the tools used in your beauty salon is essential to give customers maximum safety and hygiene. Smart Beauty Autoclave is a fully automatic and easy-to-use device with a large sterilization chamber and two removable internal trays that allow you to sterilize more tools at the same time.

Class B and compliant with European standards, this autoclave can sterilize even tools previously bagged.

The equipment is equipped with an internal printer that highlights the day and time in which sterilization is carried out, and a USB port to store the history of sterilizations made.



- ✓ Class B
- ✓ Internal printer
- ✓ USB port
- ✓ Sterilization of bagged tools



SMART BEAUTY • 45



**The autoclave water distiller is safe, reliable and easy to use. It produces 4 liters of distilled water in about 4 hours (about 20 liters of distilled water per day).**

Dimensions:  
D. 23,5 cm x W. 22 cm x H. 38 cm

### Advantages:

- Distilled water always available in the beauty center
- Use of the right water for OVERMED AUTOCLAVE
- Purchase of bulky and heavy water cans no longer needed
- Optimal maintenance of the internal components of the OVERMED AUTOCLAVE (boiler, tank, etc.)

### Technical characteristics:

- Stainless steel water tank.
- Reusable replacement filter included
- Power supply: 220-240V - 50-60hz
- Power: 750 W
- Capacity: 4 L
- Distillation volume: 1.5 L / hour
- Overheating protection: auto-power off of the distiller when reaching a temperature above 160° C.



## THERMAL BLANKET

For sauna treatments, designed to wrap the whole body, the thermo-blanket is divided into 3 independent areas (abdomen - buttocks - legs). It is equipped with an electronic controller to regulate and control the temperature, in compliance with the parameters required by current regulations, giving you the possibility to customize the treatment. A user interface, complete with LCD display, allows you to check the treatment parameters (temperature and duration) in real time. The system also provides an audible alarm of the beginning and end of treatment, as well as a warning for any malfunctions, and the display the display of total hours worked.

Available in 5 versions with different sizes.

### code 9518

Dimensions in mm open L. 1800 x W. 1700 approx.

Dimensions in mm closed L. 1900 x W. 750 approx.



### code 9533

Dimensions in mm W. 150 x L. 1500



### code 9534

Dimensions in mm W. 950 x L. 1800



### code 9520

Dimensions in mm L. 1500 x W. 500



### code 9530

Dimensions in mm L. 1800 x W. 1000





# SMART

BEAUTY



Call us on **+971 50 304 3353**

SMART BEAUTY reserves the right to make changes to the equipment described, at any time and without notice.

Smart Beauty devices are designed and programmed to obtain results on imperfections without passing the permitted aesthetic limit.

SMART BEAUTY is a Technology S.r.l. brand

